

Rechargeable Moisturization support device

MXFC-1000

Instruction Manual

Thank you for buying Maxell product.
Please read this manual to use correct in advance.
Please keep this manual carefully.

1.5hr Rechargeable

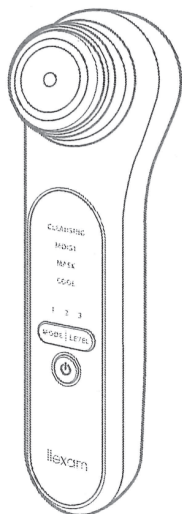





Table of contents

1	Safety Precautions	2
2	Feature of the product	5
3	Skin compatible check	5
4	Device details	6
5	Preparation before use	7
6	Example for Use	10
7	How to Use	12
8	Care of the device	22
9	How to keep	22
10	FAQs	23
11	Trouble Shooting	25
12	Specification	26
13	Disposal of product	27










1 Safety Precautions



- Precautions described here are intended for safe and correct use of this product and are to prevent injury and damage to you and others before they occur.
- Precautions are classified as follows.



Make sure to observe the precautions because they are all important in relation to safety.

 DANGER	Indicates an imminently dangerous situation in which improper handling may result in death or serious injury to a user.
 WARNING	Indicates a dangerous situation in which improper handling could result in death or serious injury to a user.
 CAUTION	Indicates a dangerous situation in which improper handling could result in slight injury to a user or damage to property (house, household effects, etc.).






<Examples of Pictorial Indication>

	<p>The  symbol indicates "Prohibited" (Do not do this).</p> <p>  "Disassembly prohibited"  "Use prohibited in wet areas"  "Use with wet hands prohibited" </p>
	<p>The  symbol indicates "Compulsory" (Make sure to do this).</p> <p>  "General action"  "Unplug" </p>





 DANGER	
 Prohibited	<p>Persons who use the following medical electronic devices should not use this product. (This product may cause the medical electronic devices to malfunction.)</p> <ul style="list-style-type: none"> - Self-contained medical electronic devices such as a pacemaker - Wearable medical electronic devices such as an electrocardiograph

 WARNING																				
 Prohibited	Do not use this product by other methods which are not described in this instruction manual. (Cause of accident, trouble, or failure)																			
	Use this product only for face and neck. (Cause of an accident or trouble)																			
	<p>The following persons should be avoid to use this product. (Cause of an accident or trouble)</p> <table border="0"> <tr> <td>- Persons with heart disease</td><td>- Persons with an acute (pain) illness</td></tr> <tr> <td>- Persons with fever</td><td>- Persons with a blood pressure disorder</td></tr> <tr> <td>- Persons with an infectious disease</td><td>- Persons with a tubercular disease</td></tr> <tr> <td>- Persons with a malignant tumor</td><td>- Persons with facial neuralgia</td></tr> <tr> <td>- Persons who are pregnant</td><td>- Persons who are menstruating</td></tr> <tr> <td>- Persons with atopic dermatitis</td><td>- Persons with an allergic predisposition</td></tr> <tr> <td>- Persons with particularly sensitive skin</td><td>- Persons with hemophilia</td></tr> <tr> <td>- Persons under dental care</td><td></td></tr> <tr> <td>- Persons whose blood vessels have expanded due to long time use of steroid hormones or liver function disorder</td><td></td></tr> <tr> <td>- Persons who are under the influence of alcohol or medication, or who are very tired</td><td></td></tr> </table>	- Persons with heart disease	- Persons with an acute (pain) illness	- Persons with fever	- Persons with a blood pressure disorder	- Persons with an infectious disease	- Persons with a tubercular disease	- Persons with a malignant tumor	- Persons with facial neuralgia	- Persons who are pregnant	- Persons who are menstruating	- Persons with atopic dermatitis	- Persons with an allergic predisposition	- Persons with particularly sensitive skin	- Persons with hemophilia	- Persons under dental care		- Persons whose blood vessels have expanded due to long time use of steroid hormones or liver function disorder		- Persons who are under the influence of alcohol or medication, or who are very tired
- Persons with heart disease	- Persons with an acute (pain) illness																			
- Persons with fever	- Persons with a blood pressure disorder																			
- Persons with an infectious disease	- Persons with a tubercular disease																			
- Persons with a malignant tumor	- Persons with facial neuralgia																			
- Persons who are pregnant	- Persons who are menstruating																			
- Persons with atopic dermatitis	- Persons with an allergic predisposition																			
- Persons with particularly sensitive skin	- Persons with hemophilia																			
- Persons under dental care																				
- Persons whose blood vessels have expanded due to long time use of steroid hormones or liver function disorder																				
- Persons who are under the influence of alcohol or medication, or who are very tired																				




⚠ WARNING

 Prohibited	<p>Do not use this product on the following regions. (Cause of accident or trouble)</p> <ul style="list-style-type: none"> - Regions where orthopedic treatment was given - Wound sites - Regions irritated due to pimples or skin rash - Regions which suffer from sensation of pain or perception disorder - Regions where metal, plastic or silicone etc., is embedded under skin - Regions where there is itching, sensation of warmth, or liver spots due to physical stimulus etc. - Regions where skin irritation is caused by cosmetics - Eyeballs and upper eyelids - Mucosal regions, such as inside the mouth - Lips
 Make sure	<p>This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.</p> <p>Children should be supervised to ensure that they do not play with the appliance.</p>
 Unplug	<p>If any abnormality occurs during charge, immediately pull out the power plug of the power adaptor from the outlet. (Cause of accident or fire)</p>
 Prohibited	<p>Use only 100-240V AC for power supply. (Cause of fire or electric shock)</p> <p>Keep away from small children. Do not use on the skin of a child. (Cause of accidents or failure)</p> <p>Never allow pins or dust to adhere to the power plug and plug for equipment. (Cause of electric shock, short circuit and ignition)</p> <p>Do not use if the power cord or the power plug is damaged or its insertion to the outlet is loose. (If the power cord breaks, a fire or electric shock may occur.)</p> <p>The power cord cannot be unnecessarily bended, pulled, twisted, bound, placed a heavy object on it, pinched, processed, or wound around the power adaptor. (If the power cord damages fire or electric shock may occur.)</p> <p>Do not use this product near inflammable substances (gasoline, benzene, thinner, etc.). (Cause of explosion or fire)</p> <p>Use only the original power adaptor. And, do not use any other product with the original power adaptor. (Possibility of electric shock or ignition)</p>
 Make sure	<p>If any unusual sound or abnormality occurs in the main body, immediately stop using this product and contact the distributor for repair together with the power adaptor. (Cause of accident or trouble)</p> <p>Follow the description of the instruction manual for the operating time and operation frequency. (Cause of accidents and trouble)</p> <p>Persons who use this product for the first time should use it while observing the condition of their skin. (Cause of skin trouble or poor physical condition)</p> <p>Use the COOL mode at an ambient temperature in the range of 15 to 35°C. (Cause of skin damage as the head is too cold if the device is used at a low temperature)</p> <p>Do not use this product beyond the specified operating time. (Cause of skin trouble)</p> <p>Do not soak the head in boiling water or dry it with hot air. (Cause of failure)</p>

⚠ WARNING

 $>\$?(\cdot)<1^*$	<p>For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance Model No. : KH-55M (Possibility of electric shock or ignition)</p>
 $C)(^*D.07^*D(0^*7\$34)^*81+7.9.0(4$	<p>Do not insert or remove the power plug in/from the outlet with wet hands. (Cause of failure or short circuit)</p>
 $C)(^*81+7.9.0(4^*3^*D(0^*\$1(\$$	<p>This mark means that this appliance should not be used near water contained in bath-tubs, basins or other vessels.</p>
	<p>Do not wash or splash water on this product. (Cause of failure or short circuit)</p>
 $B.\$)(\%9^*8^*&1+7.9.0(4$	<p>Do not use or keep this product in places easily splashed with water or in damp places, such as in a bathroom. (Cause of failure or short circuit)</p>
	<p>Do not alter disassemble or repair by unauthorized technician. (Cause of fire, electric shock, injury, or failure)</p>

⚠ CAUTION

 $-1+7.9.0(4$	<p>Check the heating dimple head, spot sensor or cooling head for any damage (scratches, etc.), and if there is any damage, do not use it. (Cause of damage to skin)</p>
	<p>Do not use this product together with other beauty treatment devices at the same time. (Cause of failure or trouble)</p>
	<p>Do not use this product if there is some injury or wound to your skin. (Cause of skin trouble or pain)</p>
	<p>Do not drop this product or give it a strong impact. (Cause of injury or failure)</p>
	<p>Do not leave this product in a place where the temperature becomes high such as the inside of a vehicle. (Cause of failure)</p>
	<p>Do not wipe with benzene, thinner or alcohol, and do not spray insecticide. (Cause of cracking, electric shock, or ignition)</p>
 $>\$?(\cdot)<1^*$	<p>Remove accessories such as earrings and necklaces, and contact lenses before use. (Cause of accident or trouble)</p>
	<p>If there is any severe skin disorder or irregularity on your skin surface, or due to downy hair, you may feel strong irritation in rare instances. In this case, avoid using this product on that region. (Cause of accident or trouble)</p>
	<p>Do not keep the cotton puff or the head on one part of your skin but always keep it in motion. (Cause of skin trouble or poor physical condition)</p>
	<p>Do not hold the cord to pull out the power plug and device plug. (Cause of disconnection, electric shock, short circuit, or ignition)</p>
	<p>When there are people around the device, pay attention to any sudden movement during use. (Cause of injury)</p>
	<p>When carrying in a bag or the like, place the product into a separate pouch. (Cause of damage to the main body or failure)</p>
 $C3\&^*<F$	<p>Make sure to pull out the power plug from the outlet while not in use. (Cause of electric shock or fire by deterioration)</p>

2 Feature of the product

- Spot sensor
Device finds spot, and notify by vibrator
- Warm function
Comfortable warm cleansing and moisturization support
- Cool function
Cool down of the skin and tighten of pores
- Memory function
Electric level is memorized after power off
- Notification timer function
Notify by "Pi,Pi" of beep sound by each approximately 45sec.
Note: No timer function on COOL mode
- Auto-off function
Leaving for a fixed time no operation, then the device shut down automatically
- Apply Titanium material on Thermal head and earthing plate
- Available to use your skin lotion and milky lotion
- Protective cap is attached

3 Skin compatible check

Please check the following to avoid your skin trouble in advance.

- 1 If you use skin lotion at first time, please check your skin by patch test.
 - (1) Please wash your arm by soap, and wipe out by towel.
 - (2) After washed your arm, please apply that new lotion to there, then keep 48hrs.
Apply a thin that lotion around coin size, and dry naturally.
Do not wet by water or scratch on that patch test place.



CAUTION

If you have abnormalities such as rash, redness, itching, irritation occur, please stop patch test immediately, and wash out.

- (3) Please confirm your skin keep normal condition after 48hrs, you can use new skin lotion.
- (4) If you have any trouble by this patch test, you can't use that new skin lotion.
- 2 Please check the compatibility between your skin lotion and this device.
 - (1) Please wash your arm by soap, and wipe out by towel.
 - (2) Please touch Thermal head around 10 sec. to your arm which is refer "STEP2 MOISTURIZE(HOT)".
 - (3) If you don't have a trouble on your skin after 48hrs, you can use that combination.

4 Device details

- ① Thermal head
(Titanium)
Vibrate while using except
"COOL" mode

- ②① Spot sensor
When Spot sensor finds a spot
or others, notify by vibration.

- ② Mode indicator

Mode	In Preparation	In Use
CLEANSING		
MOIST	Orange flashing	Orange lights
MASK		
COOL	Flashing Blue	Blue lighting

- ③ Mode switching button
Mode changes each by pressing.

- ④ Power button
Turn on by pressing the button,
and Turn off by pressing
the button again.

- ⑥ Plug socket

- ⑦ LEVEL indicator
1(Low)(Green)
2(Mid)(Green)
3(High)(Red)

After press Power button,
"1" "2" or 3 lamp lights.

- ⑧ LEVEL switching button
Electric level changes by each press.
Vibration level is no changing.

- ⑨ Charging indicator
In charging lights "Red" lamp, and
lamp for completed charging.
"Red" lamp flashing means the sigh of
charging request.

- ⑩ Buzzer
(Switching/Finishing ready/Alert/
Abnormal detective/End/
Notification timer)
Device makes Switching sound each
by pressing button, and in case of
"COOL" mode, device makes
Finishing ready sound.
In continuous using "COOL" mode is
Alert sound, and abnormal detective
sound is failure.
End sound is finishing each mode.

- ⑬ Battery(Built-in)
Lithium-ion battery : 1pc

- ⑪ Cool head(Aluminum)
Touch to the skin at
"COOL" mode

- ⑫ Earthing plate(Titanium)
Touch this place when
Thermal head is contacted
with your skin.
If your hand leaves from
earthing plate, vibration
becomes not active.

⑭

Power adaptor

Model No. : KH-55M

⑮

Power plug
Insert deeply to
the socket

⑯

**Plug to connect
the device**

⑰ Cord

⑱

Attachment

⑲ Ring(2pcs)

Ring is for cotton
fix to the device.



⑳ Protective cap

5 Preparation before use

Charging method

At first time to use after you bought this product, you need to charge the battery.


Electrical power should be used 100-240V by alternative current, then please use attached Power adaptor.

⚠ WARNING

Do not pull-out Power adaptor from socket with wet hands.
(Cause of electric shock and injuring)

Do not use in places splashed with water, such as in a bath, or in damp places.
(Cause of failures or short-circuits)


- ❶ Please turn off the power, and connect
Power adaptor to device and socket.

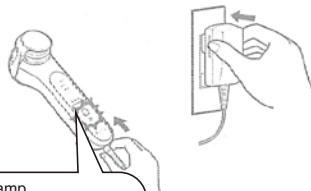

Charging lamp
Under charging : Light up red
Completed charging : Off

② After completed charging time(approximately 1.5hr), please unplug Power adaptor from the socket.

- In charging process, device and power adaptor heat up a few, but it is not abnormal.
- This device is for rechargeable exclusive. Can not use while power cord connect to the device.

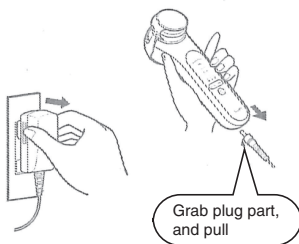

Charging time is approximately 1.5hr. approx. 7 times is available per 1 time charging.

[In case of using LEVEL 3 at all mode]

* Usage count might be reduced by circumstance temperature.

Timing of charging

Charging lamp is flashing after turn on the power button, it is the sign of charging request.

Please keep in mind the following to keep enough performance of the battery.

1. Please charge before use if you don't use 3 month or over.
2. Do not charge each by use.
3. Recommendation temperature range is 10-35deg.C. Charging time might be longer or not charge if you use this device under out of recommendation temperature.
4. Do not charge except household power source.
5. Under charging, in case of mixture noise by TV, please use another socket.

Preparation before use

⚠ CAUTION

Never use the device when you have a wound or an inflammation on your skin.
(Cause of problems with or damage to your skin.)

Remove accessories or contact lenses.
(Cause of problems with your skin or physical discomfort.)

Check the Thermal head or cooling head for any damage or deformation before use.
(Cause of damage to skin)

1 Mode and skin lotion

Mode	Cotton and skin lotion
CLEANSING(HOT) (Remove fine dirt in pores)	Please use new cotton which is soaked with your regular skin lotion.
MOISTURIZE(HOT) (Performs moisturizing care by allowing moisturizing ingredients to penetrate into the skin)	When using cotton Use new cotton soaked with your regular skin lotion. ----- When not using cotton You can also apply skin lotion directly to your skin without using cotton.
MASK(HOT) (Enhances the effect of milky lotion or sheet mask)	Use your regular milky lotion or a sheet mask by applying to your skin. Do not use cotton.
COOL (tightens the pores of your skin by cooling effect)	Apply the cooling head directly to your skin. Do not use cotton and skin lotion.


- 2** Clean skin by washing your face so that makeup and dirt will wash off.
Please wipe off your makeup with a cleansing product.

- 3** Attach cotton to the Thermal head of the device.
Remove the Protective seal from the Thermal head.
Place cotton over the head and attach it with the included ring while pressing its two projections.

- Use cotton. However, small or thick cotton may not attach to the device.

- The ring has an orientation with a clear direction.
Attach it from straight above with the wider part of the inner circumference placed toward the Thermal head.

If it is difficult to set, please attach it tightly while pressing the right and left ends (two projections) of the ring horizontally.

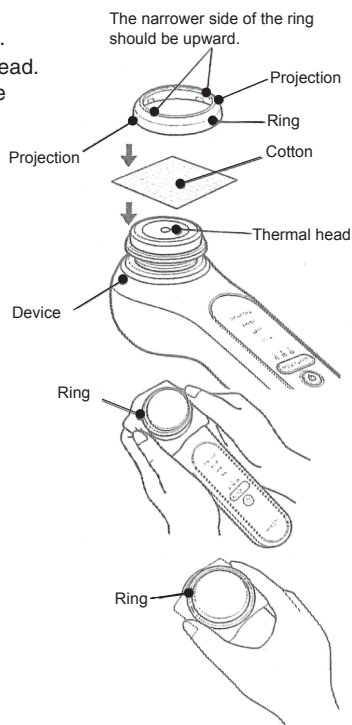


- 4** How to remove the ring.
Remove the ring by pulling it straight up while holding down the lower part.

5 Please note the following to avoid a damage your skin.

(1) Basic mode usage frequency of this device is 8 min. at once per day at a person.

- Please do not use the device for over 3 min. per day in each step. (Step 1, 2 and 3).
- Please do not continuously use the device in Step 4.

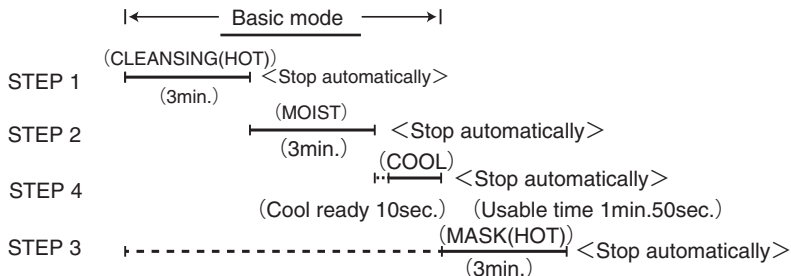

- Do not exceed to use over 11min. per day.
- At [MASK(HOT)] mode, in case of milky lotion, recommendation usage is once per day, and in case of sheet Mask, recommendation usage is once or twice per week.

(2) Please use to bare skin after remove makeup.

(3) Do not use at same place over 10sec.

(4) Please change the level to adjust your skin condition.

- You can choose "1," "2," and "3" as you like.
If your skin tingles, use the device on "1."
- If you feel something unusual with your skin, stop using the device and consult a skin care specialist.
- The device is equipped with a memory function to memorize the levels of "1," "2," and "3" that you usually use. For details, please see "About memory function" .

6 Example for Use

Press the Mode switch button and set the mode of your preferred treatment.

When using "MOISTURIZE(HOT)" without cotton, the spot sensor senses a spot.

Quick tighten-up Treatment (Required time: approximate 2 min.)

•Step 4 "COOL"

When you want to tighten your skin or cool down your skin before applying makeup, you can just use "Cool."

After 10-second preparation for cooling, it only takes one minute and 50 seconds to tighten your skin. Press the mode switch button to set it to "Cool".

Quick Treatment (Required time: approximate 3 min.)

●Step 2 "MOISTURIZE(HOT)"

When you do not have much time in the morning or when you need to go out suddenly, you can just use "MOISTURIZE(HOT)." It only takes approximate three min. to moisturize your skin thoroughly. Press the mode switch button to set to "MOISTURIZE(HOT)."

Tighten-up Treatment (Required time: approximate 5 min.)

●Step 1 "CLEANSING(HOT)"

●Step 4 "COOL"

When you want to remove fine residue remaining in pores and tighten your skin, use "Cool" after using "Hot Cleanse." This treatment tightens your skin.

After the "CLEANSING(HOT)" mode is finished, press the mode switch button to set to "COOL."

Moisturize treatment(Required time: approximate 6 min.)

●Step 1 "CLEANSING(HOT)"

●Step 2 "MOISTURIZE(HOT)"

Before you go out or when you have more time, use "MOISTURIZE(HOT)" after using "CLEANSING(HOT)." This treatment moisturizes the clean skin thoroughly.

Daily Treatment (Basic mode) (Required time: approximate 8 min.)

●Step 1 "CLEANSING(HOT)"

●Step 2 "MOISTURIZE(HOT)"

●Step 4 "COOL"

Standard daily treatment. After using "CLEANSING(HOT)" support moisture of clean skin in "MOISTURIZE(HOT)" and tighten your skin by cooling effect in "COOL"

After the "MOISTURIZE(HOT)" mode is finished, press the mode switch button to set to "COOL".

Daily Full Course (Required time: approximate 11 min.)

●Step 1 "CLEANSING(HOT)"

●Step 2 "MOISTURIZE(HOT)"

●Step 3 "MASK(HOT)" (use of milky lotion)

●Step 4 "COOL"

Care for your skin using all modes.

Fine vibrations in "MASK(HOT)" (use of milky lotion)" mode and gentle switching of the electrode will condition the skin. Tighten your skin by cooling effect in "COOL".

Weekly Full Course (Required time: approximate 11 min.)

- Step 1 "CLEANSING(HOT)"
- Step 2 "MOISTURIZE(HOT)"
- Step 3 "MASK(HOT)" (use of sheet mask)
- Step 4 "COOL"

Care for your skin once or twice a week by using all modes.

Fine vibrations in "MASK(HOT) (use of sheet mask)" mode and gentle switching of the electrode will condition the skin. Tighten your skin by cooling effect in "COOL."

7 How to Use

This device will not operate with the power adaptor plugged into an outlet.
Remove the power plug of the power adaptor from the outlet to use.

- Use after removing the protective cap.
- After the power button is pressed to "On," the Thermal head warms slowly.

⚠ WARNING

Do not use the product around eyes and ears.
(Cause of accidents or bad physical condition)

If you feel something is unusual with your skin or body while using the device, stop using it immediately and consult a skin care specialist. (Cause of accident or poor health)


Step 1 "CLEANSING(HOT)"

This mode removes fine residue remaining in pores which cannot be removed by simple washing.

- 1 Soak cotton attached to the body with skin lotion after washing your face.
Please use your regular skin lotion or rubbing-off skin lotion with an exfoliating effect.
(Do not use milky lotion or a gel-like or thick skin lotion.)
Soak the cotton with skin lotion to the extent that skin lotion permeates the back of the cotton.
Some skin lotion (just with mineral spring or water) may be difficult to bring out a good effect.

2 Press Power button

The level indicator lamp will be turned on, and "CLEANSING(HOT)" on the mode indicator lamp will start flashing. If a lamp in another mode is flashing, press the mode switch button to switch to "CLEANSING(HOT)."

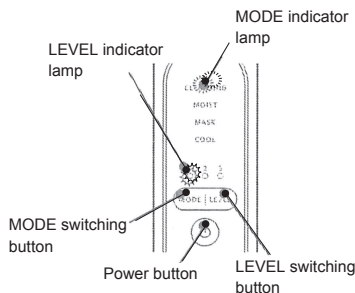

3 Select electric level by LEVEL switching button.

Switch the "strength" with the level switch button according to the condition of your skin. Every time the level indicator button is pressed, the "1", "2" and "3" of the electrical potential will be switched. There will be no change in vibration strength.

The memory function memorizes the "strength" when the power has been turned off. It is recommended that when using the device for the first time or if you feel some discomfort, use the device in "1".

4 Grab the device by hand with touching earthing plate.


Make sure to hold the body with your hand touching the earthing plate.

5 Apply cotton to the skin and move the device slowly.

Mode indicator lamp will turn on and the device will start to vibrate.

The time for use is within approximate 3 min..

Earthing plate.

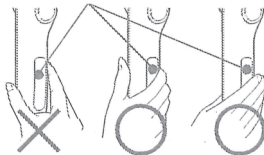

CAUTION

Unless you apply the cotton to your skin while touching the back plate, the device will not vibrate. (skin sensing function)

If the Thermal head is not vibrating, the device is not working.

If you leave the device as it is for about 3 min. without vibration, the device will automatically switch "Off." (Automatic shutoff function)

Notification timer function

Skin care time is indicated by the buzzer beeping about every 45 sec. as a guide.

[Example of skin care sequence]

Divide the area from your face to décolletage into four zones as shown in the illustration.

By moving the Thermal head in the order of ①②both cheeks->③T-zone④décolletage, the whole area can be efficiently and uniformly cared for.

Apply cotton gently to the face and slowly move the device by sliding it from the center of the face to the outer side.



⚠ CAUTION

Do not keep the cotton on the same area of skin, it is necessary to move it around. (Cause of problems with your skin or physical discomfort.)

Add more skin lotion when there is no moisture in the cotton or when the device is difficult to slide.

Do not apply the device on the same area of skin for more than 10 sec.

If you feel your skin tingling by micro current while using the device, press the level switch button and switch to "1".

If you still feel your skin tingling, use the device by adding more skin lotion.

- 6 Automatically shut down after 3min.(Auto-off function)
Notify by End sound(Pi Pi)
- 7 Remove ring and cotton
Wipe off the surface of the Thermal head with a tissue.
Do not use the cotton again once it is used. It might cause problems with your skin.

Step 2 "MOISTURIZE(HOT)"

This mode cares for skin by allowing moisturizing ingredients of skin lotion permeate into the skin.


When using cotton ※The spot sensor cannot sense the spot.

- 1 Attach new cotton to the body with the attachment ring.
For attaching a cotton puff, see "Preparation before use."
- 2 Soak the cotton with your regular skin lotion.
Do not use milky lotion or a gel-like or thick skin lotion.
Soak the cotton with lotion to the extent that lotion permeates the back of the cotton.

- 3 Press Power button.


The level indicator lamp will be turned on, and "MOISTURIZE(HOT)" on the mode indicator lamp will start flashing. If a lamp in another mode is flashing, press the mode switch button to switch to "MOISTURIZE(HOT)."

When the "CLEANSING(HOT)" mode is finished, press the power switch and the "MOISTURIZE(HOT)" mode will automatically be set.

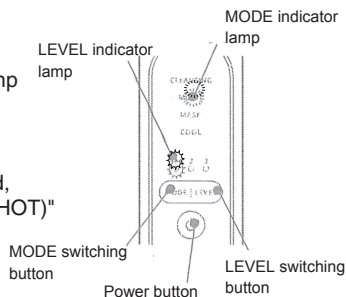

- 4 Select electric level by LEVEL switching button
Choose "the strength" of electrical potential with the Level switch button.
Switch the "strength" with the Level switch button according to the condition of your skin. Every time the level indicator button is pressed, the "1", "2" and "3" of the electrical potential will be switched. There will be no change in vibration strength.
The memory function memorizes the "strength" when the power has been shut off.
It is recommended that when using the device for the first time or if you feel some discomfort, use the device in "1".

- 5 Grab the device by hand with touching earthing plate.
Make sure to hold the body with your hand touching the earthing plate.

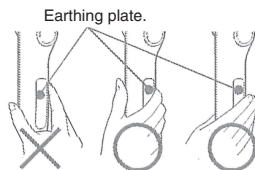

- 6 Apply cotton to the skin and move the device slowly.
Mode indicator lamp will turn on and the device will start to vibrate.
The time for use is within approximate 3 min..

CAUTION

Unless you apply the cotton to your skin while touching the back plate, the device will not vibrate. (skin sensing function)

If the Thermal head is not vibrating, the device is not working.

If you leave the device as it is for about 3 min. without vibration, the device will automatically switch "Off." (Automatic shutoff function)

Notification timer function

Skin care time is indicated by the buzzer beeping about every 45 seconds as a guide.

[Example of skin care sequence]

Divide the area from your face to décolletage into four zones as shown in the illustration. By moving the Thermal head in the order of ①②both cheeks->

③T-zone④décolletage, the whole area can be efficiently and uniformly cared for.

Apply cotton gently to the face and slowly move the device by sliding it from the center of the face to the outer side.



⚠ CAUTION

Do not keep the cotton on the same area of skin, it is necessary to move it around.
(Cause of problems with your skin or physical discomfort.)

The vibration may change while using the device. This is because the spot sensor senses the cotton. This is not a failure.

Add more skin lotion when there is no moisture in the cotton or when the device is difficult to slide.

Do not apply the device on the same area of skin for less than 10 seconds.

If you feel your skin tingling by micro current while using the device, press the level switch button and switch to "1".

If you still feel your skin tingling, use the device by adding more skin lotion.

- 7 After using 3 min., the device will turn off automatically. (Auto-off function)
Notify by End sound(Pi Pi)

- 8 Remove ring and cotton


Wipe off the surface of the Thermal head with a tissue.

Do not use the cotton again once it is used. It might cause problems with your skin.


When not using cotton

Also you can use the device without cotton in "MOISTURIZE(HOT)"

- 1 Apply a gel-like moisturizing product or thick skin lotion that you would normally use onto Thermal head and smoothen it over the entire head with your fingers.

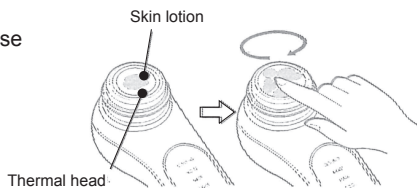

- 2 Press Power button.
LEVEL indicator will light, and "HOT MOIST" Mode indicator lamp is flashing.
In case of flashing another mode, please change MODE by press MODE switching button.
After finished "CLEANSING(HOT)" mode, it will set "HOT MOIST" mode automatically after press Power button.

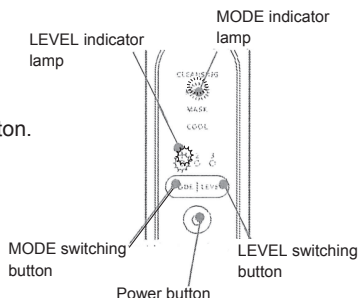


- 3 Select electric level by LEVEL switching button.

Choose "the strength" of electrical potential with the Level switch button.

Switch the "strength" with the Level switch button according to the condition of your skin. Every time the level indicator button is pressed, the "1", "2" and "3" of the electrical potential will be switched. There will be no change in vibration strength.

The memory function memorizes the "strength" when the power has been shut off. It is recommended that when using the device for the first time or if you feel some discomfort, use the device in "1".


- 4 Grab the device by hand with touching earthing plate.
Make sure to hold the body with your hand touching the earthing plate.

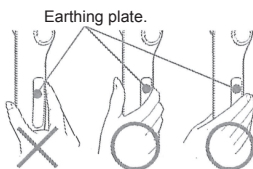

- 5 Please touch Thermal head to your skin, and move slowly.
Mode indicator lamp will turn on and the device will start to vibrate.
The time for use is within approximate 3 min.

Spot sensor

The device finds a spot on your skin, and notify by strength of vibration.

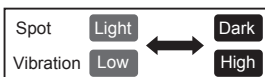

- The device notify to you by changing vibration step by step.

※ The sensor responds to dark-colored areas such as moles, shadows and bruises.

※ This device is not for Diagnosis and treatment.

Do not directly look at the spot sensor because its light is extremely bright.
Spot sensor cannot detect without moving Thermal head.

CAUTION

Unless you apply the Thermal head to your skin while touching the earthing plate, the device will not vibrate. (skin sensing function)
However, Thermal head heating start to active at he same time with turn on the power.

Notification timer function

Skin care time is indicated by the buzzer beeping about every 45 sec. as a guide.

[Example of skin care sequence]

Divide the area from your face to décolletage into four zones as shown in the illustration.

By moving the Thermal head in the order of ①②both cheeks->③T-zone④dé colletage, the whole area can be efficiently and uniformly cared for.

Apply cotton gently to the face and slowly move the device by sliding it from the center of the face to the outer side.



⚠ CAUTION

Do not keep the cotton on the same area of skin, it is necessary to move it around.
(Cause of problems with your skin or physical discomfort.)

Add more skin lotion when there is no moisture in the cotton or when the device is difficult to slide.

Do not apply the device on the same area of skin for less than 10 sec.

If you feel your skin tingling by micro current while using the device, press the level switch button and switch to "1".

If you still feel your skin tingling, use the device by adding more skin lotion.

- ❶ After using 3 min., the device will turn off automatically. (Auto-off function)
Notify by End sound(Pi Pi)

- ❷ Wipe off the surface of the Thermal head with a tissue.

STEP 3 MASK

(Regarding recommendation usage, in case of milky lotion, once per day, in case of Sheet Mask, once or twice per week is recommendation frequency)


This device provides good effective of milky lotion and sheet mask by fine vibration and slow electrode switching.

- ❶ Place a sheet mask on your face or spread milky lotion or cream that you would normally use over your entire face.

Wipe off the lotion on your hands.

- ❷ Press the power switch.

The level indicator lamp will be turned on, and "MASK(HOT)" on the mode indicator lamp will start flashing. If a lamp in another mode is flashing, please press the mode switching button to select "MASK(HOT)."
When the "MOISTURIZE(HOT)" mode is finished, press the power switch and the "MASK(HOT)" mode will automatically be set.

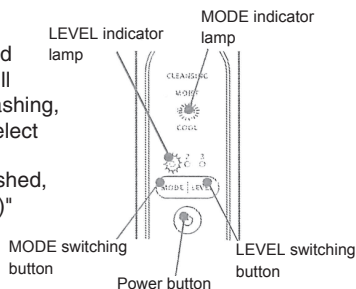


- ❸ Choose "the strength" of electrical potential with the level switch button.

Switch the "strength" with the Level switch button according to the condition of your skin. Every time the level indicator button is pressed, the "1", "2" and "3" of the electrical potential will be switched. There will be no change in vibration strength.

The memory function memorizes the "strength" when the power has been shut off. It is recommended that when using the device for the first time or if you feel some discomfort, use the device in "1."

- ❹ Grab the device by hand with touching earthing plate.
Make sure to hold the body with your hand touching the earthing plate.

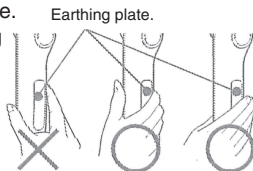

- ❺ Please touch Thermal head to your skin, and move slowly.
Mode indicator lamp will turn on and the device will start to vibrate.
The time for use is within approximate 3 min..

CAUTION

Unless you apply the Thermal head to your skin while touching the earthing plate, the device will not vibrate. (skin sensing function)
However, Thermal head heating start to active at the same time with turn on the power..

Notification timer function

Skin care time is indicated by the buzzer beeping about every 45 sec. as a guide.

[Example of skin care sequence]

Divide the area from your face to décolletage into four zones as shown in the illustration.

By moving the Thermal head in the order of ①②both cheeks->③T-zone④décolletage, the whole area can be efficiently and uniformly cared for.

Apply cotton gently to the face and slowly move the device by sliding it from the center of the face to the outer side.



⚠ CAUTION

Do not keep the cotton on the same area of skin, it is necessary to move it around.
(Cause of problems with your skin or physical discomfort.)

Do not apply the device on the same area of skin for less than 10 sec.

If you feel your skin tingling by micro current while using the device, press the level switch button and switch to "1".

- ⑥ After using 3 min., the device will turn off automatically. (Auto-off function)
Notify by End sound(Pi Pi)

- ⑦ Wipe off the surface of the Thermal head with a tissue.

STEP4 COOL

This mode provides cool feeling by the cooling head surface which is approximately 10deg.C against environment temperature, and tighten your skin by cooling effect.

⚠ WARNING

Please use 15-35deg.C temperature range.
(If you use under low temperature, cooling head will be too cool, then it is the cause of making damage to your skin)


CAUTION

Continuous use of the "COOL" mode makes it difficult to cool the cooling head.
The warning sound and the flashing "COOL" on the mode indicator lamp will inform you.
The warning sound can be canceled by pressing the power button, but the "COOL" mode cannot be used for approximately one hour after cancellation.
"CLEANSING(HOT)," "MOISTURIZE(HOT)" and "MASK(HOT)" can still be used.

"COOL" mode is high power consumption, so if battery capacity is low, it might be power off by setting "COOL" mode even other mode is available.
At that time, please charge device battery, then you will be able to use it by "COOL" mode again.

1 Press Power button.

The level indicator lamp will turn off, and "COOL" on the mode indicator lamp starts flashing. If a lamp in another mode is flashing, please press the mode switching button to select to "COOL" mode. When the " MASK(HOT) " mode is finished, press the power button, then it is set " COOL " mode automatically.

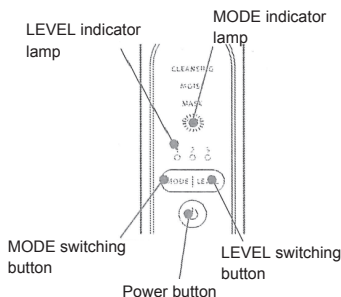

2 Wait for cool of cooling head

It is necessary 10sec. to be ready for cooling. Under preparation, "COOL" mode indicator lamp is flashing. It will be available status to use, ready end sound is made beep, then MODE indicator lamp turn on the light.

CAUTION

It is no vibration at "COOL" mode.

3 Please apply cooling head softly to your skin directly, and move slowly such as patting.

Mode indicator lamp changes to turn on the light. Cooling head does not vibrate. Available time to use is 1min. 50sec.

Since you can apply the device to the same area of your skin for only 1sec., repeatedly apply the head to your skin on and off towards the face.

Notification time function is not available on "COOL" mode.



CAUTION

Do not slide the cooling head over the skin.

Do not keep the cooling head in the same area of the skin.

The cooling effect of the cooling head may be weak.

4 After 1min. and 50 sec. past, the device will be shut down automatically.(Auto-off function)

It is notified End sound(Pi,Pi). After end of "COOL" mode, device and cooling head might be hot a few, it is not abnormal.

If remaining skin lotion/moisturizing product on the skin attached to the cooling head or device, please wipe off with a tissue.

CAUTION

If MODE indicator lamp and LEVEL indicator lamp become flashing, and the device makes abnormal detective sound ""Pi, Pi"", then it means alert of failure.

Please stop to use this device, and contact the distributor for repairing.

Skin detective function

This device is a product which operates through the human body. Applying the skin to the earthing plate and the head will make an electrical circuit start to operate.

Unless you touch the back plate and apply the head to the skin, it will not start to operate.

Memory function

This device has a memory function to memorize of the electrical strength levels of "1" "2" and "3".

In each usage mode (Step 1 to 3), the device memorize the electrical strength level just before the automatic shut down or the Power button is turned off.

If the Level switching button is pressed while the device is being used, the electrical strength level can be changed.

8 Care of the device

Wipe off any remaining cosmetic on the Thermal head, spot sensor, cooling head, main body and ring with a tissue paper, etc.

⚠ WARNING

Do not clean the main body with water.

Do not soak in water, and do not pour water on the main body. (Cause of failure)

Do not soak the head in boiling water or dry it with hot air. (Cause of failure)

⚠ CAUTION

Make sure to turn OFF the switch. (Improperly turning on the switch may give a tingling sensation or result in short-circuit.)

Do not use alcohol, benzine, or thinner. (Cause of deformation or discoloration)

Do not rub the Spot sensor with a hard brush such as a scrub brush.
(Cause of scratches and malfunctions of spot sensor)

9 How to keep

After care of the device, please attach protective cap on Warm/Cool head, avoid humidity or dust, and lay down horizontal by Thermal head upper, then keep.

Protective cap can be attached to the Warm/Cool head even with ring attached.

⚠ WARNING

Do not leave in a water splashing area or an area with significant humidity, such as bathroom. (Cause of short-circuits or failure)


Do not wind the power cord around the power adaptor or the main body.
(Cause of damage to cord)

CAUTION

Do not leave in an area subject to direct sunlight (For example, inside of a vehicle in summer) or near heaters. (Cause of deformation or discoloration)

Please attach protective cap to avoid hitting to Spot sensor.

(Cause of scratches and malfunctions of Spot sensor)

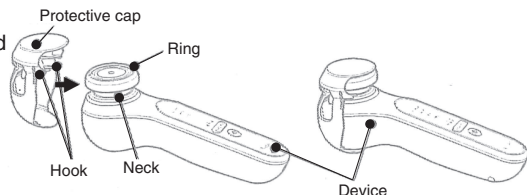

• Hook needs to adjust and attach to the neck part

10 FAQs

Please confirm the following if device is not correct movement.

If it is not returned to correct condition, please contact the distributor.

Condition	Check point	How to repair
It cannot be Power turn on.	Check Power adaptor is connected.	Please unplug Power adaptor from the device and socket.
MODE indicator lamp is not flashing.	Check the device battery capacity is empty, or not.	Please charge device battery 1.5hrs.
Do not vibrate Thermal head.	Check Power is turn on.	Please press Power button.
	Check your hand is touch to the earthing plate, or not.	It is not working Skin detection function. Please touch your hand, and apply Thermal head to your skin
	Check Skin lotion soak to the cotton, or not.	Please soak skin lotion to reach behind of cotton.
	Check Protection sheet has been removed, or not.	Please remove protection sheet.
	Check it is set "COOL" mode, or not.	It is not vibrated on "COOL" mode.
	Check your hand is touching earthing plate on "HOT CLEANSING""HOT MOIST" "HOT MASK" mode.	It is not working Skin detection function. Please touch your hand, and apply Thermal head to your skin.
Warning sound is made, and "COOL" mode indicator lamp flashing does not stop.	Check you are using "COOL" mode continuously, or not.	It is not failure. Please try to use again after 1hr where is 15-35deg.C environment. ※After press Power button and release Warning sound, "COOL" mode is not available approximately 1hr.
Power becomes off when it is set "COOL" mode	Check device battery capacity is low.	Please charge device battery.
It makes abnormal detection sound (Pi,Pi), and it does not stop all lamp flashing	Device failure.	Please contact to the distributor.
There are few available times.	Check charging condition it is out of recommendation temperature range from 10 to 35deg.C, or not.	Please charge under recommendation temperature range from 10 to 35deg.C.

Condition	Check point	How to repair
Vibration strength is not changed under Shimi-Navi detection.	Check some stain attach to Shimi-Navi sensor, or not.	Please wipe it out by tissue or cotton bar.
It cannot attach the ring.	Size of cotton is different each by kind, then it might not be suitable to the ring size.	Please use other cotton.
	The ring is set upside down.	Please attach narrow side of inside up to cover the cotton.
Cooling head does not become cool.	Check continues using by "COOL" mode.	Please try to use it again after 1hr at 15-35deg.C environment.
	Check environment temperature is over 35deg.C.	You might feel cool that is different between cooling head surface and skin temperature. Please use it at 15-35deg.C environment temperature.
Warm head vibrates even that does not touch to your skin.	Check you touch the earthing such as water supply.	It might be vibrated if you touch the earthing plate through the earth. Please do not touch the earth.
Warm head is not heating.	Check environment temperature is less than 15deg.C.	When environment temperature is low, heating up time might take a time. Please use it at 15-35 deg.C environment temperature.
	Check environment temperature is over 35deg.C.	You might feel warm that is different between Warm head surface and skin temperature. Please use it at 15-35deg.C environment temperature.
It cannot charge to the device battery.	Check charging condition it is out of recommendation temperature range from 10 to 35deg.C, or not.	Please charge under recommendation temperature range from 10 to 35deg.C.
The device charging time is long.	Check charging condition it is out of recommendation temperature range from 10 to 35deg.C, or not.	Please charge under recommendation temperature range from 10 to 35deg.C.
The device turns off when you are using.	Check the device and Power adaptor connect, or not.	Please unplug Power adaptor from the device and socket, then use it.
	Check the device battery capacity is empty.	Please charge device battery 1.5hrs.

11 Trouble Shooting

Q1. Tingling sensation is felt on the skin.

- Check the intensity level. If it is set on "3" level, please press the level switching button to change to "1" level.
- If skin lotion in the cotton pad wears off, replenish it again.
- ※ If you still feel uncomfortable on the skin, please stop using the product immediately and consult a dermatologist.

Q2. Can "CLEANSING(HOT)" mode remove makeup?

- This product cannot remove makeup.
Please use a makeup remover cleanser.

Q3. Will dirt be more easily removed effective by applying more stronger force?

- Electric power is used for removing dirt and moisturizing. Thus, by applying more stronger force, it will not make any difference.
Gently apply to skin and massage slowly.

Q4. Cotton pad cannot be attached properly.

- Use a normal rectangular or round cosmetic cotton pad.
(If the cotton pad is too thick or thin, it will not be able to attach.)

Q5. Can the product be used in the bathroom?

- This product is not waterproof.
Never use or place this product in wet and damp areas, e.g. bathroom.

Q6. Does the "COOL" mode work with skin lotion?

- Yes, but it will not enhance the cooling effects.

Q7. The cooling head does not feel cool during the "COOL" mode.

- The "COOL" mode automatically stops after being used around 2 min.
Approximately 10 sec is required at the start of the mode before the cooling effect takes place.
It will last for another 1 min 50 sec before it stops automatically.
- The "COOL" mode may be used continuously and continuous use will reduce the cooling effect.
Do not continuously use the "COOL" mode.
When the cooling head is continuously applied on the same area of skin, the temperature of the skin surface decreases and the cooling effect of the head may be lessened. Check it by applying the cooling head on another area.

Q8. Can I use milky lotion or gel?

- If you use cotton, you can not use cosmetics other than skin lotion.
Please use skin lotion which has keratin removal effect or others.
- Milk-based lotion or gel can be used when using " MOISTURIZE(HOT)" mode without cotton pad and "MASK(HOT)" mode.
Please refer to the "About mode switching" in the "Preparation before use" section in the instruction manual for details.

Q9. Does this product have whitening or spot removal function?

- This product is not for whitening or removing blemish.
It will pull out fine dirt that remains deep in the pores and penetrate moisture into the skin.

Q10. Which regions are not suitable to use this product?

- Please refer to the instruction manual for the list of regions.

Q11. Another case

- If you still feel uncomfortable on the skin, please stop using the product immediately and consult a dermatologist.
- If the device has any abnormal condition, please contact to the distributor which is you bought.

12 Specification

Rated voltage	100 - 240V AC, common to 50-60Hz
Rated Power Input	9W (Power Adaptor : KH-55M)
Power consumption	Approximately 8W (Charged AC100V AC)
Battery	Li-ion battery (3.7V DC 1620mAh 1pc)
Charging time	1.5hr (Room temperature 10-35deg.C)
Rated time	Approximately max.11 min. (Basic mode : approximately 8min.)
Weight	Device : Approximately 210g (not including Protective cap and Ring) Power adaptor : Approximately 75g
Attachments	Ring (2 pcs), Protective cap (1pc)

13 Disposal of product

DANGER

- With the removed batteries, do not charge and avoid the following.
(Cause of heat generation, ignition, or explosion)
Do not throw batteries into a fire, or overheat them.
Do not puncture with a nail, give an impact, disassemble, or modify them.
Do not bring (+) and (-) into contact with metal, etc.
Do not carry or keep them together with necklaces, hairpins, etc.
Do not leave batteries in a high-temperature area, such as near a fire or in a burning sun.

WARNING

- Keep the removed battery out of the reach of children.
If it is suspected that a child has swallowed a battery, seek for medical consultation immediately.
- If the battery fluid leaks, do not touch it with bare hands but take the following measures:
If the fluid gets in your eyes, do not rub eyes but immediately wash them with clean water, and then seek for medical consultation.
If the fluid gets on your body or clothing, thoroughly wash it away with clean water, and then seek for medical consultation.
- Do not disassemble except when disposing of the product.
(Cause of fires, electric shocks or injury)

The battery is to be disposed of safely.

■ This product uses a lithium ion battery.


Lithium ion battery is recyclable, valuable resources. When disposing of the product, handle the lithium ion battery by following your local regulation.

In addition to the battery, other recyclable parts are included. When disposing of parts other than the battery, follow your local regulation.


How to remove Lithium ion battery

- ① Hold the power switch and level switch buttons simultaneously for 5 sec. or more, and all indicator lamps light and then the mode and level indicator lamps light sequentially.
Leave it until they turn off.


- 2 Remove the three screw concealing pieces and the three screws on the back of the body.

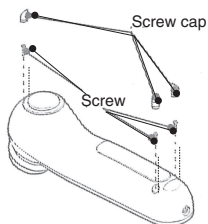

- 3 Insert a flat blade screwdriver etc., into the gap between the body cases (front) and (rear), and remove the body case (rear).

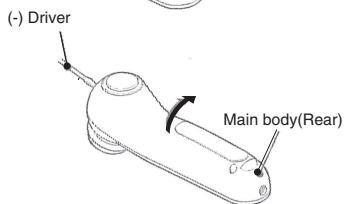

- 4 Disconnect the lithium ion battery connector, and take it out of the battery holder.

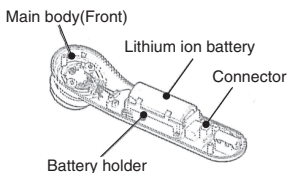


⚠ DANGER

Do not connect (+) and (-) terminal of removed Lithium ion battery.
(Cause of short circuit and burn of skin)

CAUTION


Make sure to pull out the device plug from the main body after charging.
Do not store to keep wrapping around by cord to the body of device and also avoid to keep connecting device plug.

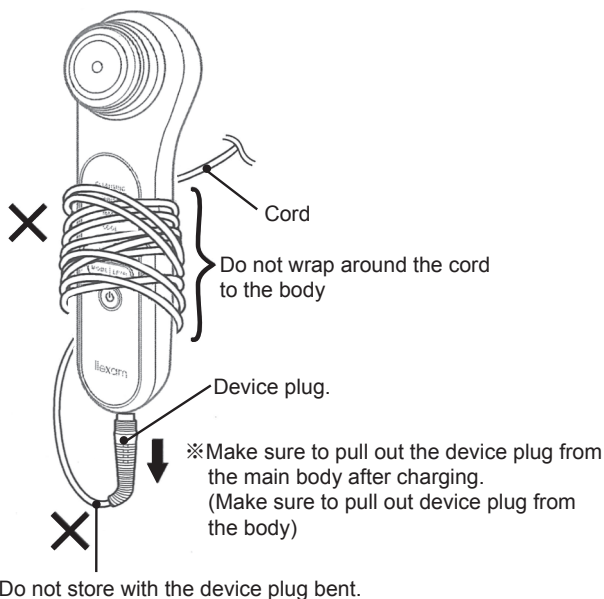

NOTICE

Attention for Use

- Make sure to hold the main body of the product with your hand touching back plate and apply the heating dimple head of the product to your skin. Otherwise, the product will not be operated.

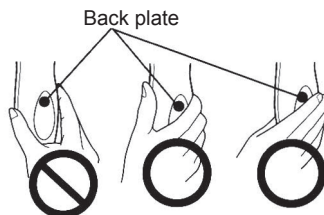

About Skin Detective Function

- This product which is operated through the human body.
- The electric circuit operates on skin in contact with the back plate and the heating dimple head of the product.
- The product may keep on operating even after the heating dimple head is moved away from skin, but this is not a failure. In such a case, it can be stopped by turning off the power switch or releasing the back plate.

About Operating Time

- The continuous operating time of each mode is default.
 - The time is not included under below circumstances:
 - When the heating dimple head or the head attached with cotton puff and back plate are not contacted with skin continuously.
 - The cotton puff does not contain enough toner.
- * **The actual operating time normally would be longer.**

About Auto-off Function

- If the product is not in operation around 3 minutes, the power will turn off automatically.
- There is an ending sound given and the product will turn off automatically when each mode is completed.

Maxell, Ltd.

30-13 Motoyoyogi-Cho, Shibuya-Ku, Tokyo 151-8527

<http://www.maxell.co.jp>